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User Manual

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1	Functional elements	26
2	General	27
2.1	Foreword	27
2.2	Scope of delivery	27
2.3	CE mark	27
2.4	Manufacturer's declaration	28
3	Introduction	29
3.1	General information about blood pressure	29
3.2	Intended Use	29
4	Operating the device	30
4.1	Special features of the PressureTel™	30
4.2	Put the batteries into the PressureTel™	30
4.3	Setting time and date	31
4.4	Installing BodyTel™ Mobile software	31
4.5	Set up your BodyTel™ user account	33
4.6	Pairing the PressureTel™	34
4.7	Putting on the PressureTel™	35
4.8	Measuring blood pressure and saving the values	36

4.9	Calling up previous blood pressure values	36
4.10	Transmitting values	37
4.11	Deleting stored the values	37
4.12	Legend for symbols	38
5	Important instructions	39
5.1	Safety measures and important instructions	39
6	Troubleshooting	40
6.1	Error messages	40
6.2	Error messages of the PressureTel™	41
6.3	Troubleshooting	42
7	Care and maintenance	43
7.1	Cleaning	43
7.2	Metrological inspection and adjustment	43
8	Warranty	44
8.1	Warranty	44
8.2	Repair	44
9	Technical data	45
10	Manufacturer's electromagnetic guidelines	46

1 Functional elements

- ① "START/STOP" button referred to as "left button" in the text
- ② "MENU" button referred to as "right button" in the text
- ③ Bluetooth® interface
- ④ Cuff with cuff tube
- ⑤ Tube connection
- ⑥ Display
- ⑦ Cuff tube

The battery compartment is to be found at the bottom side.
(see chapter 4.2)

2 General

2.1 Foreword

The integrated autofeedback system allows you to perform first-class, comfortable and extremely gentle blood pressure measurements.

Measured values can be transmitted automatically via a compatible mobile phone* to your personal online patient file of the *BodyTel Center* (www.bodytel.com).

Thank you for your trust. We are available for service and product questions.

*Infos: www.bodytel.com, support@bodytel.com or free BodyTel™ Support +800-0-2639835.

2.2 Scope of delivery

1. PressureTel™ blood pressure monitor for home use with cuff (size M) and cuff tube
2. Operating instructions
3. Storage bag
4. Two batteries
5. Warranty card

2.3 CE mark

PressureTel™ complies with the provisions of Council Directive 93/42/EEC of 14 June 1993 as well as national laws concerning medical devices.

PressureTel™ complies with the provisions of Council Directive 99/5/EC (R&TTE). It carries the CE mark as proof of this. PressureTel™ complies with the provisions of the European Standard EN 1060-1 Non-invasive sphygmomanometers - Part 1: General requirements, and EN 1060-3: Supplementary requirements for electro-mechanical blood pressure measuring systems.

2.4 Manufacturer's declaration

Manufacturer: I.E.M. GmbH
Cockerillstraße 69
52222 Stolberg, Germany

I.E.M. GmbH declares that the product type designation: "Stabilograph mobil®" is in accordance with the relevant provisions of Council Directive 93/42/EEC of 14 June 1993 concerning medical devices and Council Directive 1999/5/EC of 9 March 1999 concerning telecommunications terminal equipment.

Stolberg, August 2010, **I.E.M. GmbH**

Attention: BodyTel Europe GmbH markets the product **StabiloGraph™ - mobil** under the name of **PressureTel™**.

Measuring accuracy

The measuring accuracy is in accordance with the regulations of the European Society of Hypertension (ESH) and the British Hypertension Society (BHS).

Compatibility and costs

PressureTel™ can only be used with **compatible mobile phones**. A list containing those phones can be found on www.bodytel.com.

To always transmit measured values and manually entered events to your patient file, you either can use a WLAN or you need to use the **mobile internet**. In that case please make sure you have access to this service and check your mobile phone plan for related costs. Maybe consider a data plan.

3 Introduction

3.1 General information about blood pressure

- Your blood pressure fluctuates throughout the day. A normal maximum for adults would be a blood pressure of 140/90 mmHg.

(Recommendations for measuring blood pressure, German Society of Hypertension, Heidelberg, 2001)

- If you monitor your blood pressure regularly, you will discover more about your body's reactions to external influences (e.g. stress or exercise) and will be able to develop a healthier lifestyle.
- Please take your PressureTel™ with you to each doctor's appointment. In case your doctor does not have an internet connected PC in his office, you may bring along a printed version of your log book. Your doctor will be able

to evaluate your data and will thereby be in a better position to care for your health.

3.2 Intended Use

- PressureTel™ is intended for home monitoring of blood pressure and pulse.
- Never change your medication on the basis of results you have obtained. The evaluation of the measurements and the treatment based on them may only be carried out by the attending physician.

4 Operating the device

4.1 Special features of the PressureTel™


Please once carry out steps 4.2 and 4.3 below before starting. In case the device has already been preconfigured (batteries are already inserted and time and date is already set), please go straight to point 4.4.

i Caution: In order to immediately transmit measured values to your online logbook you need have your mobile phone's Bluetooth® interface activated, the *BodyTel Mobile* application running and must be close to the PressureTel™ (max. 10 meters).

In case you need help, please contact the free BodyTel™ Support +800-0-2639835.

4.2 Put the batteries into PressureTel™

Open the battery compartment at the bottom side of the casing. Insert the batteries supplied with the device, as shown below. Make sure that the ribbon is underneath the batteries!

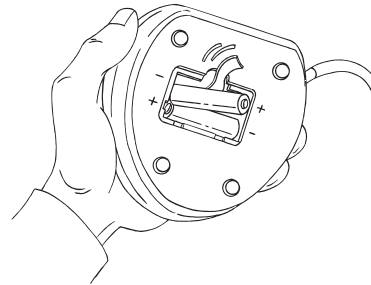4 Operating the device

4.3 Setting the time and date

Hold the **right button (MENU)** until MON, DAY, HOUR and MIN are flashing in the display. Next, press the **left button (START/STOP)** repeatedly until the correct minute (MIN) is set. Press the **right button** again to call up the hour (HOUR) and then the date (YEAR, MON and DAY). Use the **left button** to set the correct value in each case. The setting ends when no button has been pressed for 30 seconds or by pressing the **right button** until the year no longer appears.

ⓘ Attention: Please adjust the device's time only, if all measurements taken on the device have already been sent to the *BodyTel Center* (via mobile phone, PC or Gateway).

4.4 Installing BodyTel™ Mobile

1. Activate *Bluetooth®* on your mobile phone and be sure the Bluetooth® visibility is set to "**visible to all**".

2. You may choose between two options to download the BodyTel Mobile application to your mobile phone:

Option A - Open the internet browser on your mobile phone and enter **<http://mobile.bodytel.com>**.

- OR -

Option B - Use your home computer, enter **www.bodytel.com** in the address field of your internet browser, go to "**Service & Downloads**" and choose "**Downloads**". Now enter your mobile phone number in the respective field. BodyTel™ now sends you **a text message to your mobile phone** that includes the download link. Open the text message on your phone, navigate to the link and open it. **This will start your phone's internet browser.** Please accept all connection requests that might come up.

3. As soon as the BodyTel™ download page appears on your mobile phone screen, please choose your language by clicking the respective link.
4. Please read and accept BodyTel™'s „**Privacy Statement**“ as well as the „**Business Terms**“ by checkmarking the two respective fields and click the „**send**“ button.
5. Please choose now the respective software version according to your mobile phone and confirm and accept all following download and installation requests.
6. When asked please choose the phone's **internal memory** as the installation location.
7. The *BodyTel Mobile* software will now be downloaded and installed.
8. Please start now the *BodyTel Mobile* application on your phone. Several connection and access requests might possibly show up. Please always allow maximum access by choosing "**Always allowed**".

① **IMPORTANT:** Some mobile phones save the rights you grant right away, other phones need to be configured manually in order to stop asking over and over again. In order to assign the rights manually please follow these instructions:

- Open the phone's menu and navigate to the item "**App. Mgr.**".
- In this application manager please navigate to "**BodyTel**" (depending on your phone you might need to choose "**Installed Programs**" first).
- Click "**open**" or - depending on your phone - click Options and choose "**Settings**" or "**Suite-Settings**" to get to the correct menu.
- Assign "**Always Allowed**" to the following points: "**Network Access**", "**Connectivity**", "**Read user data**" and "**Edit user data**".

4 Operating the device

- Additionally, some phones also allow to set a fixed „**Access Point**“ as well as the behavior of the *BodyTel Mobile* „**Auto-start**“ in that menu.

For help please check your mobile phone user manual or call the BodyTel™ Customer Support at +800-0-2639835.

4.5 Set up your BodyTel™ User Account

1. Open the *BodyTel Mobile* application and the Installation Assistant welcomes you. It will guide you through the process so please just follow the instructions on your mobile phone's screen.
You are going to be asked for a user name (a valid email address), a personal password (at least 6 characters long) and some other data.

2. In the end you will be sent an **activation email** to the entered email address that contains an **activation link**. Please check this email with your home PC und click on the activation link provided in that email. This will automatically open your PC's webbrowser. Please now follow the instructions on your personal BodyTel™ login site in order to complete your profile.

This is mandatory in order to be able to use all BodyTel™ system functions.

3. You can now turn your attention back to your mobile phone and your PressureTel™.

4.6 Pairing the PressureTel™

1. Please start the *BodyTel Mobile* application on your mobile phone.

❗ **Attention:** If your mobile phone again asks you several connection and access requests, please close the application and follow the instructions of the grey box at the end of chapter 4.4.

2. When starting the *BodyTel Mobile* application for the first time, a message will show up, asking, if you would like to pair a device now. Please confirm with "**Yes**".

❗ **Attention:** In case you accidentally chose "**No**" or you want to pair a device later you can also get to that menu by choosing "**Options**" - "**Settings**" - "**Device Manager**".

3. Please now choose the **PressureTel™** from the device list and follow the instructions (**4. - 11.**).

❗ **Attention:** Your mobile phone should be very close to the PressureTel™ because the PressureTel™ will always find the counterpart with the strongest *Bluetooth®* signal. Make sure that no other *Bluetooth®* enabled devices (other mobile phones, computers, etc.) are close to the PressureTel™ while pairing.

4. Place your PressureTel™ next to your mobile phone and hold the **right button** (MENU) of the PressureTel™ until the display changes.
5. Now press the **right button** of the PressureTel™ again several times until you see "**P02**" flashing in the Pressure-Tel™ display.
6. Press the **left button** (START/STOP) of the PressureTel™ and please pay attention to your mobile phone.
7. Watch your mobile phone's display and confirm and accept all pairing requests that might come up on it.

4 Operating the device

8. Now your phone will ask for a *Bluetooth®* password.
9. Please immediately enter the password **6624** *on your mobile phone* and click "**OK**". A single short "beep" *on the PressureTel™* will confirm the correct code.
10. After the "beep" please immediately perform a blood pressure measurement as described in chapters 4.7 and 4.8.


① **Attention:** This measurement is part of the pairing process and will NOT appear in your logbook.

11. After another single "beep" of the PressureTel™ the *BodyTel Mobile* application on your phone will state that the pairing process is now finished.

① **Note:** A successfully paired device will be listed on the *Device List* in the *BodyTel Mobile* application.

4.7 Putting on the PressureTel™

Assemble the upper arm cuff as shown, with the label on the outside. Now pull the cuff over your left arm, until the lower cuff edge lies approx. 2 cm above the elbow joint. The artery marking should be on the inside of your elbow. To produce correct values, you should keep your arm, hand and fingers still during the measurement.

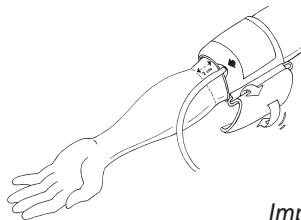
Important note:

We recommend carrying out blood pressure measurements after approx. 3–5 minutes sitting at rest.

4 Operating the device

4.8 Measuring blood pressure and saving values

Wearing the cuff as described in chapter 4.7 press the **left button** (START/STOP) briefly to start the measurement. The PressureTel™ inflates and then deflates in stages. After successful measurement the display will show (example):

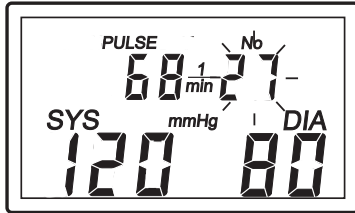The reading number (here 27) flashes. The measurement is stored automatically. If you do not want to store the value, press the left button within 10 seconds.

4.9 Calling up previous blood pressure values

First, hold the **right button** (MENU) until the display flashes. You will see the average value of the last 3 measurements, shown by a flashing "A".

Then press the **left button** (START/STOP). The reading and the date/time of the last-stored measurement will now be shown alternately.

Press the **left button** repeatedly to see the last (max. 99) values on the PressureTel™.

Press the **right button** briefly or wait 30 seconds without pressing a button, if you do not wish to call up any further blood pressure values.

This is only possible in case you already have measured values in your PressureTel™. If not the display will immediately go to the edit Time and Date state after you pressed the **right button**.

4 Operating the device

4.10 Transmitting values

A data connection sends the value measured to your on-line logbook at www.bodytel.com via your mobile phone. If previous measurements have not been transmitted, all unsent values will be transmitted together.

When data are being sent via Bluetooth®, "**C02**" will be shown on your PressureTel™. After a successful transmission the blood pressure monitor will emit a brief acoustic signal (1 x beep). After each blood pressure measurement the value is transmitted automatically.

If necessary, values can also be transmitted manually. Therefore, press the **right button** (MENU) of the monitor (first time: long), until "**C02**" flashes. Then press the **left button** (START/STOP) button once. "**C02**" stops flashing and the stored values are transmitted.

ⓘ Attention: To transmit data from the PressureTel™ to your mobile phone, the *BodyTel Mobile* application must be opened.

4.11 Deleting stored values

If you press the **left button** (START/STOP) for longer than 5 seconds an acoustic signal will be emitted and the display will show "**Clr**".

If you press the **left button** for further 5 seconds, all previously stored values will be deleted and the time mode will reappear.

4.12 Key to symbols

Key	Operation / situation	Measure / solution
♥	Pulse detected	
SYS (mm/Hg)	Displays systolic value (upper value)	
DIA (mm/Hg)	Displays diastolic value (lower value)	
PULSE (1/min)	Displays pulse (Pulse rate per minute)	
HOUR	Hour	
MIN	Minute	
DAY	Day	
MON	Month	
	Batteries are empty	Replace batteries (see 4.2)
Err 01	Error message	See 6.1
No	Value number	
Clr	Delete values	If the left button is pressed for further 5 seconds, all measurements are deleted
PO2	Pairing	See 4.6
CO2	Mobile data transmission	See 4.10

5 Important instructions

5.1 Safety measures and important instructions

Please take note!

If you experience discomfort during a measurement – press the **left button** (START/STOP) to deflate the cuff immediately. Open the cuff and take it off. Please contact BodyTel Europe GmbH.

In addition, please note:

- If you suffer from heartbeat irregularities, you should consult your doctor about whether blood pressure monitoring with the PressureTel™ is appropriate for you.
- Please also consult your doctor if you are taking anticoagulants (e.g. Marcumar, Aspirin)
- To produce precise values, the device should stand on a hard surface. The cuff should be fixed around the bare upper arm (see diagram). For measuring, sit comfortably on a chair with your hands resting on your thighs.
- The device is suitable for upper arms with a circumference of 24–32 cm when using the M-cuff and 32–38 cm when using the L-cuff. The L-cuff can be obtained directly from BodyTel Europe GmbH.
- Do not under any circumstances use the device in the bath or the shower. If, however, liquid gets into the device, you must remove the batteries immediately and refrain from further use. Contact our service department.
- The device may not be used near a nuclear magnetic resonance tomograph.
- Using components that are not part of the scope of delivery can produce faulty measurements (Exception: battery replacement).

6.1 Error Messages

Error code	Cause	Solution
Err 01	Movement artefacts	Keep your arm still and repeat the measurement See 6.2, 4.3, 4.7
	Pressure remaining in cuff	Repeat measurement only after 10 seconds
Err 02	Cuff	a) Put on the cuff correctly (see 4.7) b) Upper arm circumference 24–32 cm for the M-cuff and 32–38 cm for the L-cuff
	Different patient	c) Start new measurement
Err 03	System error	a) During measurement: repeat the measurement
	Pneumatic fault	b) After replacing the batteries: reinsert batteries, contact the service department
Err 04	Battery fault	a) Replace the batteries (see 4.2)
		b) Rechargeable batteries used (Very limited number of measurements)

6 Troubleshooting

6.2 Error messages PressureTel™

Message	Cause	Solution
Cod 01	a) Pairing error b) Mobile phone switched off	Try pairing again (see 4.6) Switch on mobile phone
Cod 02	a) Connection problem b) No values in the memory c) All values already transmitted	a) Manual transmission (see 4.10) b) Performing a measurement c) Performing a measurement
Cod 03	Meas. could not be sent	manual transmission (see 4.10)
Cod 04	Communication error	if the problem recurs please contact the manufacturer
Cod 05	Communication error	if the problem recurs please contact the manufacturer
Cod 06	Hardware communication error	if the problem recurs please contact the manufacturer
Cod 07	Memory error	if the problem recurs please contact the manufacturer
Cod 08	Undefined error	if the problem recurs please contact the manufacturer

6.3 Troubleshooting

Description of error

During the measurement process, pressing the left button (START/STOP) once does not stop the measurement immediately and the cuff does not immediately deflate	Open the cuff, remove the batteries and do not use the device further. Contact the service department.
Left button pressed, but the desired measurement does not take place (pumping stops prematurely)	Check the batteries and replace them if necessary
No symbols appear in the display	Check batteries and replace them if necessary
Battery symbol appears in the display	Replace batteries
Time and/or date are incorrect	see 4.3: Setting the time and date
Faulty measurements	a) Put on the cuff correctly (see 4.7) b) Upper arm circumference of 24–32 cm for the M-cuff and 32–38 cm for the L-cuff c) Keep your arm still during the measurement

7 Care and maintenance

7.1 Cleaning

Cleaning: If required, wipe the plastic casing and the cuff with a damp cotton cloth only.

Important note: The PressureTel™ cuff may not be washed or dry-cleaned.

7.2 Metrological inspection and adjustment

The manufacturer has calibrated the device for a period of two years. Metrological inspection and adjustment must take place at least every two years. Metrological inspection and adjustment is subject to a fee and can be carried out by I.E.M. GmbH, a responsible authority, or an authorised maintenance service.

Disposal of used electrical and electronic appliances, including batteries:

This symbol on a product or its packaging indicates that this product must not be disposed of with normal household waste, but must be handed over to a designated collection point for the recycling of electrical waste and electronic equipment. For more information, please contact your local city office, your household waste disposal service or the shop where you purchased the product.

8.1 Warranty

- 1.** BodyTel Europe GmbH grants the blood pressure monitor a warranty of two years from date of sale. Proof of purchase must be supplied in the form of a correctly completed warranty card or a receipt.
- 2.** Defects due to faulty material or workmanship will be repaired free of charge during the warranty period.
- 3.** The warranty period is not extended as a result of warranty services being provided, either for the device as a whole or for any components.
- 4.** The warranty excludes:
 - a)** any damage caused by improper use, e.g. by failure to observe the operating instructions.
 - b)** damage caused by maintenance or interventions carried out by the purchaser or unauthorised third parties.
 - c)** damage resulting from transport, either from the manufacturer to the customer, or on the way to customer services.

d) component parts subject to normal wear and tear (cuff, batteries etc.)

- 5.** Liability is excluded for direct or indirect damage caused by the device, even if the damage to the device falls within the terms of the warranty.
- 6.** Unauthorised modification of the warranty card invalidates it.
- 7.** Further claims, for whatever reason, are excluded.

BodyTel Europe GmbH does not warranty the batteries supplied.

8.2 Repair

In the case of a warranted repair, please complete the warranty card, giving precise details of the sender and the defect. Please contact our service department (see front cover), who will inform you about shipping possibilities.

9 Technical data

Measuring range	Systole (SYS): 70 to 260 mmHg Diastole (DIA): 30 to 180 mmHg Pulse (PULSE): 40 to 240 (1/ min)	Data transmission	Bluetooth Class II
Operating temperature	+10 °C to +40 °C (+ 50 ° F to +104 ° F)	Power supply	2 alkaline 1.5 V AA batteries (AA, LR6)
Operating humidity	15% – 90%	Unit dimensions	125 mm x 131 mm x 51 mm (W x D x H) 4,9 " x 5,2 " x 2,0 " (W x L x H)
Storage temperature	-20°C bis +50°C (-4°F to +122°F)	Weight	270 g (9,52 oz) incl. batteries
Storage humidity	15% – 95%	Technical modifications reserved. (DEU)	
Indicating range	Pressure: 0 to 300 mmHg Pulse: 40 to 240 beats per minute (1/ min)		
Accuracy	Pressure: ± 3 mmHg in indicating range		
Measuring method	Oscillometric		
Cuff size	Upper arm circumference: M-cuff 24–32 cm (9,5 " to 12,6 ") L-cuff 32–38cm (12.6" to 14.96")		
Memory capacity	99 measurements		
Display	Liquid crystal display (LCD)		
Operating elements	Two buttons: »START/STOP« and »MENU«		

10 Manufacturer's electromagnetic guidelines

Guidelines and manufacturer's declaration – electromagnetic immunity

The PressureTel™ blood pressure monitors are intended for use in the electromagnetic environment specified below. The customer or the user of the PressureTel™ blood pressure monitor should ensure that it is used in such an environment.

Emissions test	Compliance level	Electromagnetic environment – guidelines
RF emissions CISPR 11	Group 1	The PressureTel™ blood pressure monitors use RF energy only for their internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The PressureTel™ blood pressure monitors are suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.

Recommended safe distances between portable and mobile high-frequency telecommunication devices (e.g. mobile phones) and the PressureTel™ blood pressure monitors

The PressureTel™ blood pressure monitors are intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the PressureTel™ blood pressure monitor can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the PressureTel™ blood pressure monitor, according to the maximum output power of the communications equipment, as shown below.

Safe distance as a function of transmission frequency in m			
Rated output of HF device in W	80 MHz to 800 MHz	800 MHz to 2.5 GHz	150 KHz to 80 MHz
	$d=1.17 (P/W)$	$d=2.33 (P/W)$	
0,01	0.12	0.23	not applicable
0,1	0.37	0.74	not applicable
1	1.17	2.33	not applicable
10	3.70	7.37	not applicable
100	11.70	23.30	not applicable

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

10 Manufacturer's electromagnetic guidelines

Guidelines and manufacturer's declaration – electromagnetic immunity

The PressureTel™ blood pressure monitors are intended for use in the electromagnetic environment specified below. The customer or the user of the PressureTel™ blood pressure monitor should ensure that it is used in such an environment.

Immunity test	IEC 60601-1 test level	Compliance level	Electromagnetic environment – guidelines
Electrostatic discharge (ESD) IEC 61000-4-2	+/- 6 kV contact discharge +/- 8 kV discharge in air	+/- 6 kV contact discharge +/- 8 kV discharge in air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial and hospital environment.

10 Manufacturer's electromagnetic guidelines

Guidelines and manufacturer's declaration – electromagnetic immunity

The PressureTel™ blood pressure monitors are intended for use in the electromagnetic environment specified below. The customer or the user of the PressureTel™ blood pressure monitor should ensure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment –guidelines
			Portable and mobile communications equipment should be used no closer to any part of the PressureTel™ blood pressure monitors including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance: $d > 0.3 \text{ m}$
Radiated RF interference IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	$d = 1.17 (P/W)$ for 80 to 800 MHz
			$d = 2.33 (P/W)$ for 800 MHz to 2.5 GHz
			Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment that is marked with the IEC symbol for non-ionizing radiation.

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Field strengths from fixed transmitters, such as base stations for radio (cellular) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the PressureTel™ blood pressure monitors are used exceeds the applicable RF compliance level above, the PressureTel™ blood pressure monitor should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the PressureTel™ blood pressure monitor.

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Important: In case of warranty please fill out all sections of this card and send it with the device.

Rev. 27.03.2012
BPT03000021

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Copyright © Juli 2010 BodyTel Europe GmbH, alle Rechte vorbehalten.